

2000

Examination
and procedures **chair**




Rehance®

2000

Examination and Procedures Chair

The Model 2000 is a multi-purpose chair that may be used in examination and screening rooms, or for laser procedures.

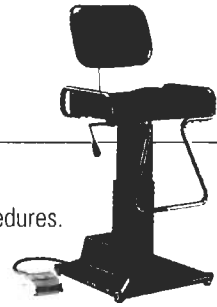

Power Base

A corded foot switch, which raises and lowers the Model 2000, may be placed at a convenient location around the chair's base. Seat height ranges from 19 3/4" fully lowered to 28 3/4" in the fully elevated position.

Designed for Comfort

The Model 2000 is equipped with handrails and a footrest for patient security and comfort. The 10" x 15" upholstered back adjusts both vertically and horizontally. An 18" square base provides exceptional stability.

Rotation Lock

The Model 2000 swivels 360°. The chair top is locked in place with a rotation lock lever on the right-hand side of the seat.

Upholstery

Black upholstery is standard, but you can match your office or treatment room décor with our selection of 14 easy-to-maintain vinyl colors. Other vinyl colors and nylon tweed inserts are available at additional cost. All materials are fire retardant. To request a color selection card, contact your dealer or call Reliance at 1-800-735-0357.

Dimensions and Specifications

Base: 18" square. Seat height: 19 3/4" lowered, 28 3/4" raised. Approximate shipping weight: 98 lbs. Electrical: 115 volts, 60 Hz, 2.4 amperes.



Precision by Tradition