




# Endurance™

Instrument Stand

*Unparalleled strength & style.  
Complement your Endurance™ Chair.*

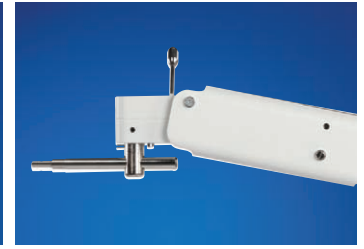

The Reichert® **Endurance™ Instrument Stand** meets the needs of the eye care professional who demands the very best in durability, functionality and style.

## Endurance™ Instrument Stand Features:

- **Strong, slim and stylish.** Durable construction outlasts the competition while contemporary styling complements your exam room. The slim profile and small footprint allows easy access to your patients and instruments without compromising floor space.
- **Logical control panel.** Buttons are located in close proximity to each other for easy access to the main on/off power switch, chair elevation and control, lower instrument arm, overhead lamp, projector and accessory instruments.
- **Convenience for accessories.** Three rechargeable instrument wells with charge indicators. Corded (BIO) instrument power and voltage selector. Three additional accessory outlets allow for consolidation of power cables.
- **Functionality at its finest.** The Phoroptor® arm comes fully counterbalanced for a Reichert® Ultramatic RX™ Phoroptor® and the lower instrument arm is counterbalanced for the Reichert® XCEL® 255 Slit Lamp. The built-in counterbalance adjustment mechanisms can be easily accessed in each arm. The height range on the lower arm accommodates wheelchairs and standard patient access.



Reichert® **Endurance™** Tilt Chair & Stand

## Endurance™ Instrument Stand Specifications

**Catalog Number** 15253

### Dimensions

Height to top of pole: 68.5 in. (174 cm)

Base: 21.75 in. W x 22 in. D (52 cm x 56 cm)

### Weight

240 lbs. (109 kg)

### Refractor Arm

Range in vertical motion: 12.5 in. (32 cm)

Length fully extended\*: 26.25 in. (64 cm)

\*Including Phoroptor® pole mount

Counterbalance weight limit: 20 lbs. (9 kg)

### Lower Arm

Range in vertical motion: 24 in. to 37.5 in.

(61 cm x 95 cm)

Length fully extended: 35 in. (89 cm)

Counterbalance weight limit: 60 lbs. (27 kg)

### Electrical

120 VAC, 50 - 60 Hz

### Fuses

6 Amps Main

### Corded Instrument Voltages

6 V, 7.5 V, 12 V, 12 V Fixed

## Optional Accessories

- 15254** Counterbalanced Keratometer Arm
- 15255** Non-counterbalanced Keratometer Arm
- 15258** Projector Mount
- 15002** Gas Lift Stool
- 15003** Gas Lift Stool with Backrest